
Two cameras simultaneously recording

- https://github.com/yangll0620/multiCams/releases/tag/v1.0
- Installation
 - Download the source code (zip) and unzip it by right click -> Extract All...



- Usage
 - Launch the Anaconda Prompt through Start menu -> Anaconda3(64bit)-> Anaconda Prompt (anaconda3)
 - ➡ If no, install anaconda here https://docs.anaconda.com/anaconda/install/windows/
 - In the terminal
 - Open the multiCams folder by typing cd yourmultiCamsFolder and Enter
 - Activate multiCams environment by typing 'conda activate multiCams' and Enter
 - Start program by typing 'python multiCams.py -n 2', 2 represents simultaneously record 2 webcams
 - Anaconda Prompt (anaconda3)
 - (base) C:\Users\lingling>cd C:\Users\lingling\Desktop\multiCams-1.0\multiCams-1.0
 - (base) C:\Users\lingling\Desktop\multiCams-1.0\multiCams-1.0>conda activate multiCams
 - (multiCams) C:\Users\lingling\Desktop\multiCams-1.0\multiCams-1.0>python multiCams.py -n 2_

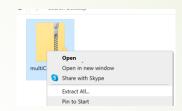